


Cryotherapy - laryngology

Cryo-S Electric II

- Nose bleeds
- Trigeminal neuralgia (TGN)
- Tonsillitis chronica
- Nasal polyps
- Snoring
- Nasal congestion
- Other diseases of the mucous membranes of the upper respiratory tract and skin
- Cryoablation of detritus stasis



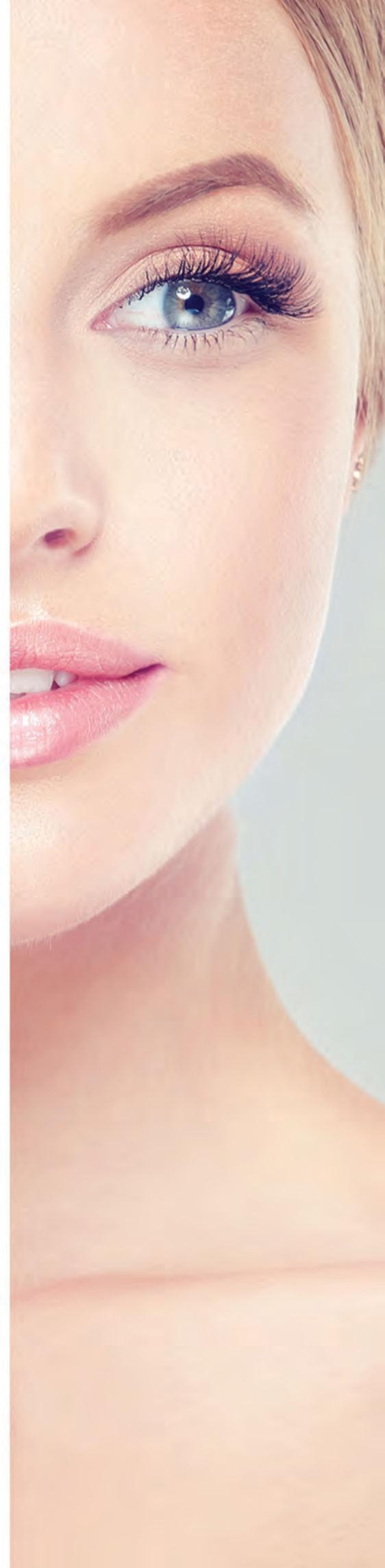


Enlarged nasal concha blocking nose



Nasal concha after cryotherapy



Appearance of the tonsils and throat before cryotherapy

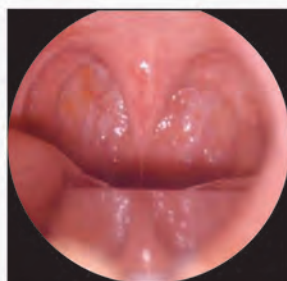

Condition after cryotherapy

Cryotherapy ENT

Cryo-S Electric II - state-of-the-art cryosurgical device manufactured by METRUM CRYOFLEX is the new generation of apparatus generally approved by many experts in the field from 1992. Cryo-S Electric device is controlled by micro-processor and all the parameters can be monitored on a LCD touch screen.

Cryo-S Electric II smart technology - the new standard in cryosurgery. What is new?

- Mode selection, probe cleaning and freezing can be performed automatically using foot switch or touch panel so the site of a lesion remains sterile.
- Electronic communication (chip system) between the unit and all kinds of probes connected to it. Unit recognizes optimum operating parameters and auto-configures to probe characteristics.
- The gas and flow pressure is set automatically, any manual adjustment is not necessary
- Contains automatic probe cleaning system and two freezing modes (continuous and discontinuous)
- Temperature of the probe, pressure, gas flow, time of the procedure are screened during freezing.

CRYOTHERAPY ADVANTAGES:

- Office procedure
- No scar formation
- No primary bleeding (pre & post cryo)
- Palliation of painful neoplasms
- Well accepted by elderly & debilitated patients
- No unpleasant odour of burning tissue (As in SMD)
- Absence of scarring as associated with deep cautery
- Post cryo recovery is smooth
- No major complications encountered
- Most cost effective as compared to laser or radiofrequency therapy

MAIN INDICATIONS:

- Nasal bleeding (epistaxis)
- Trigeminal neuralgia (tgn)
- Tonsillitis chronica
- Nasal polyps
- Snoring
- Nasal congestion
- Other diseases of the mucous membranes of the upper respiratory tract and skin
- Cryoablation of detritus stasis



Contact probes for Cryo-S Electric II. All of the probes are single use or reusable, easy to sterilize and easy in exploitation.

DISTRIBUTOR

BVM MEDICAL LTD
 BVM House, Trinity Lane, Hinckley, LE10 0BL,
 United Kingdom
 Tel: 01455 614 555 Web: www.BVMMedical.com

PRODUCER

METRUM CRYOFLEX
 ul. Kolejowa 16A
 05-092 Łomianki
 POLAND, EU

export@metrum.com.pl
www.metrum.com.pl

